

making life

easy

# Transit Wheelchair

User Manual



# User Manual

Read all operating and safety instructions before operating your wheelchair.

## Assembly Instructions

### To Unfold:



**01**  
Remove the chair from the box and place it on the ground. Remove any cardboard protection.



**02**  
Lift folding push handles up so that the release pin engages the upper arm.



**03**  
Place both hands on the upholstery just to the inside of the seat tubes and gently push down and out until the seat is in a horizontal position.



**04**  
**BE SURE BOTH HANDS AND ALL FINGERS ARE FLAT ON THE UPHOLSTERY.**  
Do not obstruct the tubes that secure the seat. Continue pushing downward and outward until the tubes fully engage with the seat supports.

WARNING: FAILURE TO FOLLOW STEP 4 COULD LEAD TO POSSIBLE HAND OR FINGER INJURY.

### To Fold:




**01**  
Pull upward on the middle of the seat from the front and rear until the chair is completely folded.



**02**  
Release the folding push handles by pushing down on the plastic release paddle and allowing the push handles to fall into a folded position.

### To Attach Legrests:

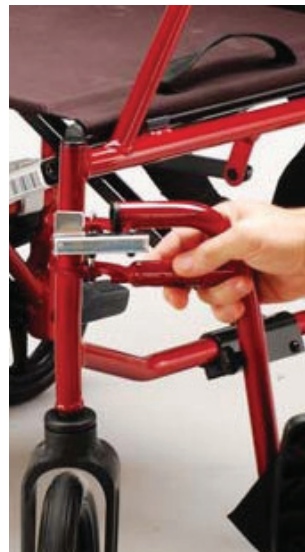

**01**  
First determine the left and right side by folding down the footplate and placing them in front of the seat. The footplate's should be facing each other.

**02**  
Place the holes above the pegs and slowly lower the holes onto both pegs.

**03**  
Swivel the Legrests inward until the silver release paddle engages and locks the footrest into the forward position.

### To Remove The Legrests:

**01**  
Push the silver release paddle inward to release the footrest from the forward position.

**02**  
Swivel outward and lift the Legrests from the footrest pegs.



## Safety and Maintenance

Wheelchair must be in a stationary position, on a level surface with the brakes locked when sitting or attempting to sit. Risk of serious damage to the wheelchair and serious injury may result if brakes are not in the locked position.

Do not exceed the maximum user weight. This may cause damage to the wheelchair, and injury to the user, as the wheelchair can become unstable and tip over.


This wheelchair was designed to carry one person at any time. Carrying more than one person at any time may exceed the weight limit and make the chair unstable and therefore unsafe.

Any damage caused to the wheelchair due to miss use or exceeding the maximum weight limits is not covered by warranty.

This wheelchair was designed to be easy to use on even surfaces. Should you come across surfaces with uneven, bumpy or rough ground we advise to avoid if possible.

Even though this wheelchair is ideal for indoor or outdoor use, we recommend that you avoid outdoor use during bad weather (such as rain or snow) as this may increase the risk of incidents. Exposure to excess moisture or water may cause rust to the steel frame.

We recommend that the maximum incline angle to climb with your wheelchair is 10%. Any attempts to climb or descend steeper inclines may cause the wheelchair to become unstable and therefore unsafe.

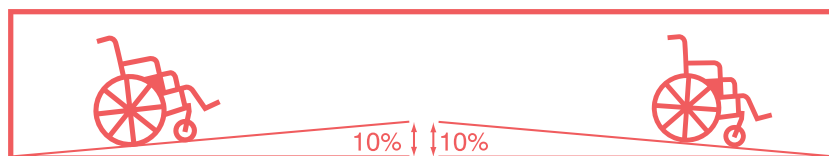

### Maximum Safe Incline Angle (Ascending and Descending)

Curbs and other obstacles should be approached with caution. It is recommended that the wheelchair is driven around such obstacles wherever possible. If these obstacles cannot be avoided we recommend these are approached slowly and at 90° from the front or rear of the wheelchair.



**Correct Curb Approach**



**Incorrect Kerb Approach**

It is very important to regularly check and service the wheelchair, ensure the brakes are working correctly and that all of the nuts and bolts are secure. Do not use the wheelchair if the brakes are not working correctly. You should keep the wheelchair clean and free from dirt, as this will prevent dirt and grim from causing damage to the wheel bearings, castor bearings and hinges.

Worn tyres may make your wheelchair difficult to drive. Ensure that both the rear drive wheel and front castors are not excessively worn. If the wheels are excessively worn, dented, twisted or out of alignment we recommend that these are replaced.

Please contact your place of purchase to arrange the purchase of any spare parts you need.

## Specification

Colour Burgundy	Overall Width 630 mm
Seat Width 460 mm	Gross Weight 11.6 kg
Seat Depth 400 mm	Net Weight 10.3kg
Seat Height 480 mm	Folded Dimensions 255x675x760 mm
Back Height 440 mm	Max User Weight 100 kg
Arm Height 210 mm	Warranty 12 months
Push Handle Height 910 mm	

## Product Warranty

### Goods with warranty term

- Goods have a 12 month warranty from the date of purchase.
- Valid proof of purchase is required for all warranty claims.
- The warranty does not extend to any subsequent purchaser or owner of the goods.
- The warranty covers manufacturing faults only, and does not cover change of mind or damage, fault or failure due to external causes, including (but not limited to) accident, abuse, misuse, incorrect assembly or normal wear and tear.
- Any unauthorised repairs, alterations or modifications or any attempt to open or dismantle the good(s) will render the warranty void.
- Defective goods, together with valid proof of purchase, should be returned to the place of purchase or shipped directly to the manufacturer. If shipped directly to the manufacturer, you must prepay the shipping and transportation charges and insure the shipment or accept the risk of loss or damage during shipment and transportation.
- The manufacturer will determine at no cost to you the cause of the fault (if any) within a reasonable time frame. Any repaired or replaced goods, original goods determined not faulty, or refund (if a major failure is determined) will be sent to the pharmacy you returned the goods to for collection, or shipped to you freight prepaid if you shipped the goods directly to the manufacturer.
- Our goods come with guarantees that cannot be excluded under Australian Consumer Law. You are entitled to a replacement or refund for a major failure and you are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure. Goods presented for repair may be replaced by refurbished goods of the same type rather than being repaired.



making life  
**easy**

---

Ph: 1300 363 303  
[www.makinglifeeasy.com.au](http://www.makinglifeeasy.com.au)  
Level 8 637 Flinders Street  
Docklands VIC 3008